

Non Clip Clasby Humeral Brace

Material - Bioprene



Manufactured
in the UK




Designed in
Association with
Clinicians



Hand Wash &
Air Dry Only



Same Day
Dispatch

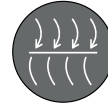

Breathable



ISO:13485
Accredited

Manufactured in the UK. Clinically proven, the non clip clasby is effective at controlling pain and is associated with high patient satisfaction scores. Innovative clam shell design gives a universal fit for the left or right. Chest strap clip sewn in one end to provide comfort and ease when tensioning strap. Constructed from Bioprene which is a new generation of anti bacterial material treated with AEGIS microbe shield. It provides compression for full contact support. Also available with plastic stays.



Indications for use

- » Mid and proximal diaphyseal fractures

Features and benefits

- » Maintains integrity of the skin more than traditional rigid braces
- » The chest strap distributes weight across the body whilst not impinging on the breast area, especially in females
- » It supplies traction and provides mild abduction to control and manage pain
- » Easy to apply, enhancing compliance and promote early rehabilitation
- » There are 6 easy to remove metal stays that exert pressure equally to stabilise the fracture
- » Trimmable around the elbow to provide a more customised fit
- » There are 2 dual pull compression straps to apply compression and prevent rotation of the fracture
- » Extension piece provided to accommodate oedema.



Sizing Information

Size	AC Joint Elbow	Code
Small	300	CBNC/3
Medium	320	CBNC/5
Large	340	CBNC/7
X - large	350	CBNC/8