

UNIVERSAL Non Clip Clasby Humeral Brace

Material - Bioprene



Manufactured
in the UK




Designed in
Association with
Clinicians



Hand Wash &
Air Dry Only



Same Day
Dispatch

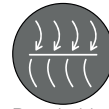

Breathable



ISO:13485
Accredited

Manufactured in the UK. Clinically proven, the Universal non clip clasby is effective at controlling pain and is associated with high patient satisfaction scores. Innovative clam shell design gives a universal fit for the left or right. Chest strap clip sewn in one end to provide comfort and ease when tensioning strap. Constructed from Bioprene which is a new generation of anti bacterial material treated with AEGIS microbe shield. It provides compression for full contact support.



Indications for use

- » Mid and proximal diaphyseal fractures

Features and benefits

- » Maintains integrity of the skin more than traditional rigid braces
- » The chest strap distributes weight across the body whilst not impinging on the breast area, especially in females
- » It supplies traction and provides mild abduction to control and manage pain
- » Easy to apply, enhancing compliance and promote early rehabilitation
- » There are 6 easy to remove metal stays that exert pressure equally to stabilise the fracture
- » Trimmable around the elbow to provide a more customised fit
- » There are 2 dual pull compression straps to apply compression and prevent rotation of the fracture
- » Extension piece provided to accommodate oedema.



Sizing Information

Size	AC Joint Elbow	Code
Small	300	CBNC/3
Medium	320	CBNC/5
Large	340	CBNC/7
X - large	350	CBNC/8